

Portable Massager


INSTRUCTION MANUAL



SPECIAL FEATURES

- No peeling, painless, no side-effects.
- Relieves fatigue, relaxation and promotes blood circulation.
- It is convenient - Massage any part of the body.
- Adjustable massage speed.
- It's modern figure with clean lines and the combination of traditional and contemporary design boundaries have created a product of aesthetic excellence.
- Advanced technology ensures safe operation.
- Independent design of the circuit board provides double safety protection.
- It's easy-to-operate feature enables the smooth operation. It is the best choice for wellness centre's, beauty salons, massage parlours, slimming centre's and home use.

COMPONENTS

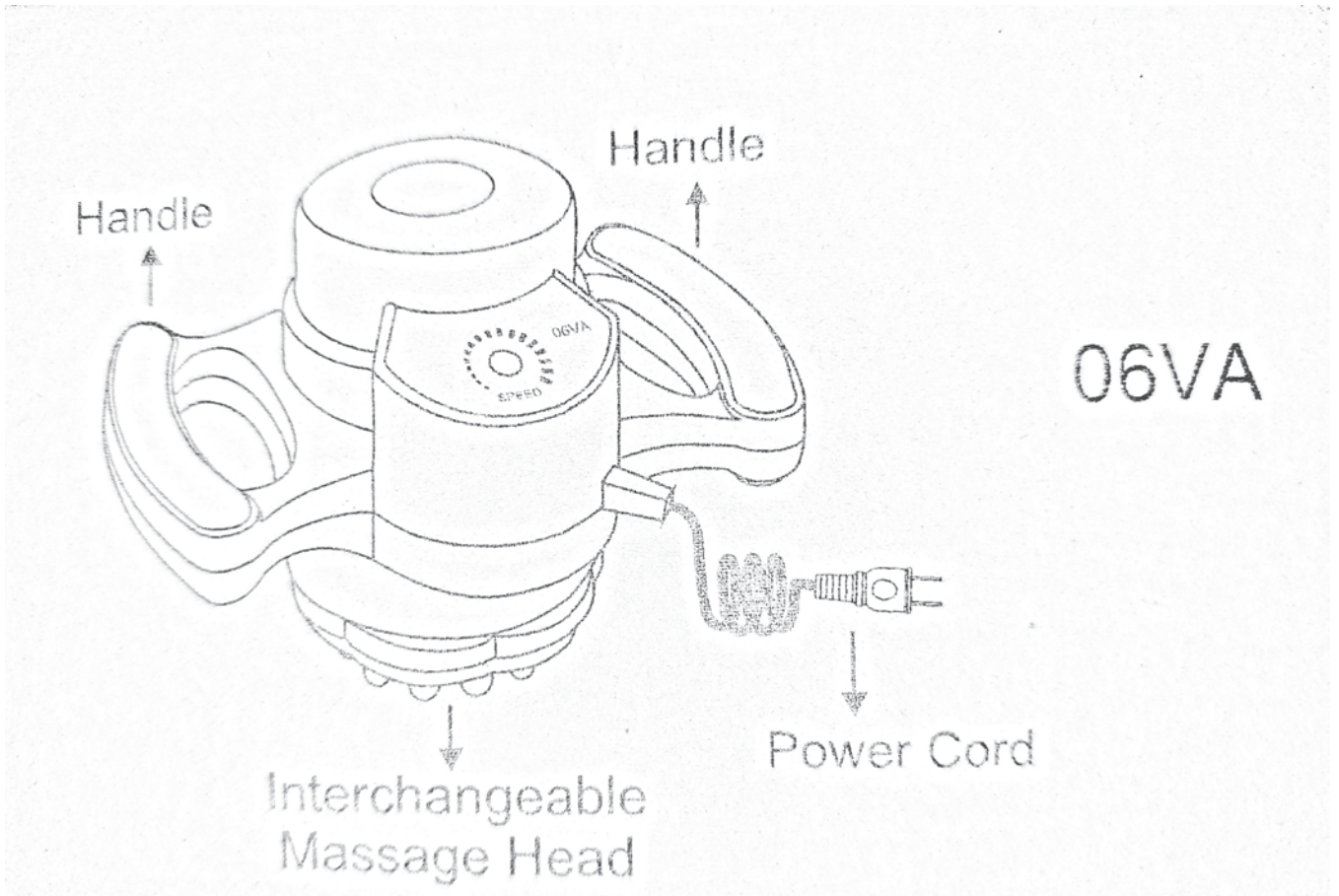

Interchangeable Massage Heads:



Circle-dots Massage Head




Half-sphere Massage Head



Columnar Massage Head

OPERATING INSTRUCTIONS

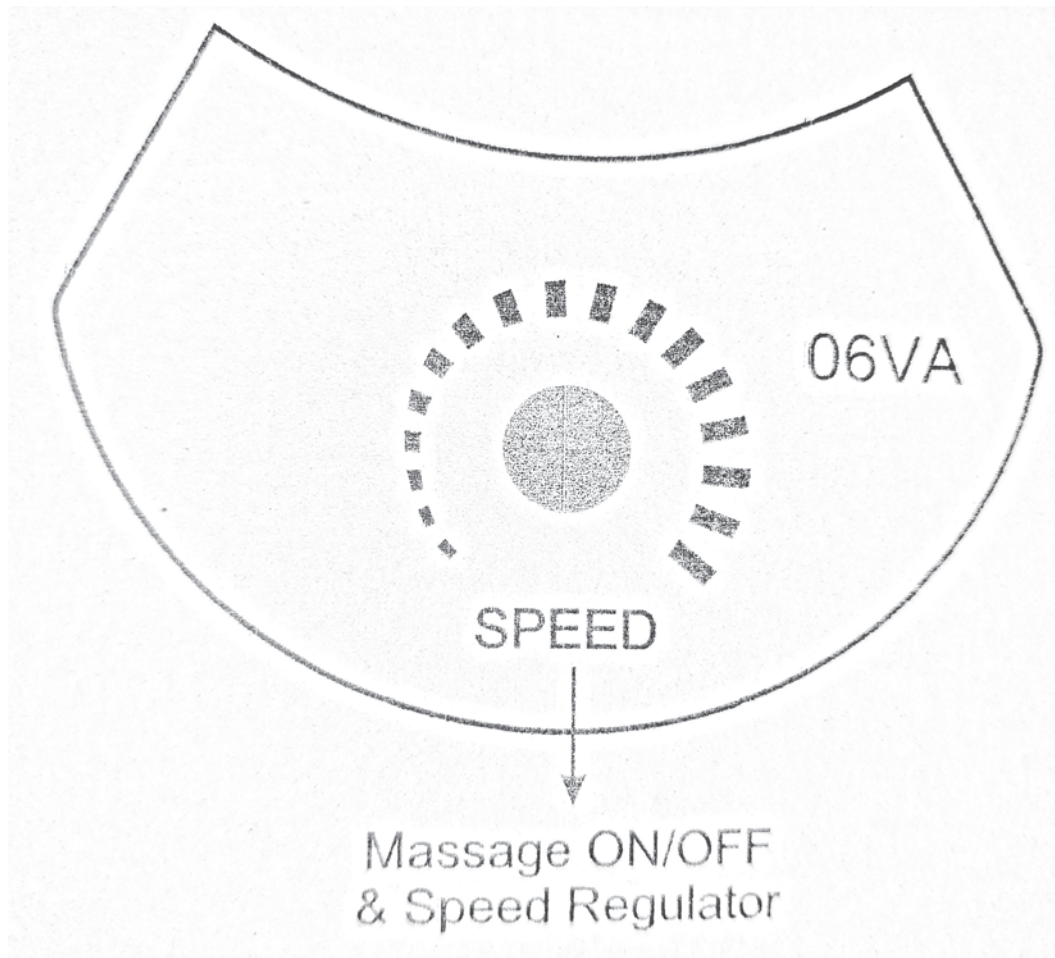

- 1. Connect the unit to the electrical supply.**
- 2. Assemble the desired massage head.**
 - Three interchangeable massage heads attachments are available to be interchanged to the unit according to desire.
- 3. The treatment area should be covered with cotton cloth instead of synthetic fiber which may cause skin sensitivity under high speed vibration.**
- 4. When treatment with gel or oil, please wrap the massage head to prevent dirt from getting into the internal section, resulting in damage to the unit.**
- 5. Hold the handles and turn on the speed regulator in accordance to the desired.**

SAFETY INSTRUCTIONS

1. Do not start the unit if it is not ready, OR it may increase rubbing abrasion of the massage head.
2. Clean the unit and accessories with clean water or light detergent instead of the chemical substances.
3. When in the process of treatment, make sure to avoid injured skin.
4. Person undergoing the following diseases should NOT use this unit:
 - Children
 - Pregnant
 - Women Menstruating
 - Acute inflammation
 - Thrombophlebitis
 - Lymph disease
 - Lung Hypoxia,etc.
5. If the supply cord is damaged, it must be replaced by the manufacturer, service agent or similarly qualified persons in order to avoid a hazard.
6. The unit is driven by motor. The brush inside the motor will be rubbing-abrasion by high speed running. So the brush needs to be replaced with a new one after normal working for about 600 hours. If in overload working for a long period of time, it may need to be replaced earlier. The brush to be replaced should be the original parts and it should be processed by authorised technician.

Contact Global Wellness Enterprises if any problems occur.

SPECIFICATIONS

| | |
|-------------------------|--------------------------|
| Product No. | 06 VA |
| Rated Voltage | 220V |
| Rated Frequency | 50Hz |
| Rated Power Consumption | 100W |
| Weight | 2.5kg |
| Package Dimensions | 28.3cm x 25.3cm x 20.5cm |